


# *Dental* **MadaJet XL**

*Needle-Free Injector for Patient Comfort!*

*The  
New Generation  
Because Patients  
Hate Needles*

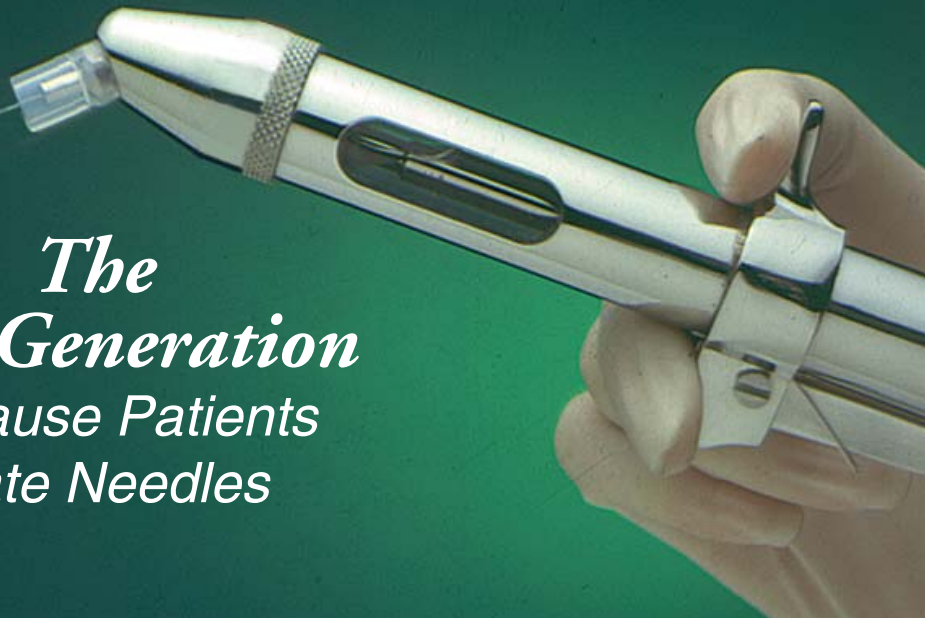

*The MadaJet XL is a proven means  
of reducing fear and trauma associated with  
traditional needle and syringe injections:*



- *Injects local anesthetics into palatal and soft tissues with greater ease and virtually no pain.*
- *Provides instant local anesthesia or regional blocking for many dental and Maxillofacial procedures.*
- *Excellent patient acceptance for over 25 years.*
- *Has interchangeable Extenda Tips for easy sterilization.*
- *Can be used with all age groups.*
- *Comes with a full year warranty.*
- *Provides adequate anesthesia for deep needle insertion.*



**MADA, INC.**  
DOMESTIC/INTERNATIONAL  
625 Washington Avenue  
Carlstadt, NJ 07072  
Tel: 201-460-0454 • Fax: 201-460-3509

To Order Call Toll Free: **800-526-6370**

## Specifications:

Consistent depth of penetration of 2-2.5mm below epithelium and makes a wheal at the base of injection of 5-6mm in diameter (Not Topical).

Consistent volume of 0.1cc per injection intradermally.

Permits approximately 38 injections with single loading (to 4.0cc).

May be sterilized by autoclaving or your usual sterilization process -  
**DO NOT USE DRY HEAT.**



### **Package No. 301 Contains:**

one MadaJet XL, two 1/2" (127mm) Extenda Tips®, holder, wrench, stylets, two extra pyrex fill-chambers, custom carrying case and MadaCide-FD Disinfectant/Cleaner. Additional Extenda Tips® available.

1 1/2" Extenda Tips® available on request.

**MadaJet XL Dental  
Technical Presentation  
available FREE upon request.**

## Applications:

- Extraction or preparation of deciduous teeth.
- Curettage and scaling.
- Mental and nasal palatine blocks.
- Adequate anesthesia prior to deep-needle insertion for palatine and infraorbital blocks.
- Adequate infiltration for long buccal blocks, without need for further anesthesia.
- Mental foramen blocks via infiltration at buccal surface of mandible.
- Minor gingivectomies.
- Incision and drainage.
- Pulp anesthetic infiltration.
- Pointing of abscesses.
- Copper tube impressions.

### **Maxillary Procedures**

- Similarly for posterior superior blocks for upper molars, middle superior blocks of bicuspids and posterior palatine blocks.

### **Mandibular Blocks**

- Two injections prior to deep needle insertion for mandibular blocks will eliminate pain from the needle and fluid pressure.