



# Mada Oxygen Regulators



R1835GBHF-25

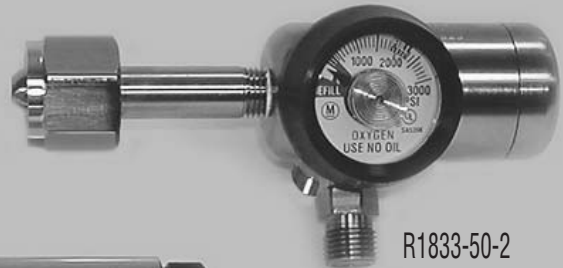

R1833-50-2

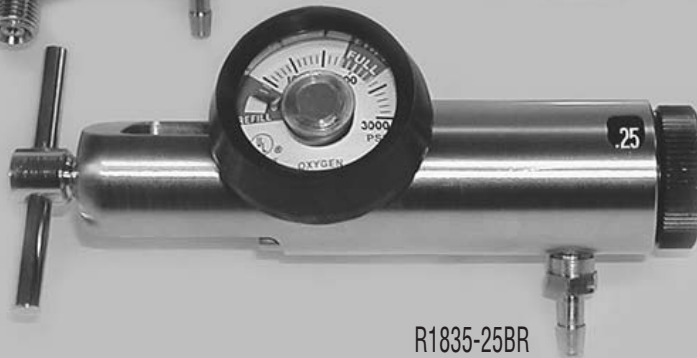

R1835-25BR

*Mada...Quality You Can Count On!*

*Mada offers a wide range of oxygen regulators that are designed for both hospital and EMS applications.*

*Mada oxygen regulators offer the features you expect, including easy to read flow rates and a dial-in click-style flow controller at the most competitive prices.*

*Best of all, Mada regulators offer you a five-year manufacturer's warranty that you can depend on.*

## **Features**

- 100% Brass construction for durability and safety.
- Large, easy to read flow rate numbers.
- Compact, low-profile contents gauge with gauge protector.
- 5 Year Warranty that you can count on.
- Large variety of regulators available, as well as accessories.



**MADA, INC.**

DOMESTIC/INTERNATIONAL  
625 Washington Avenue, Carlstadt, NJ 07072  
Tel: 201-460-0454 • Fax: 201-460-3509

To Order Call Toll Free: **800-526-6370**

Web Page: [www.madamedical.com](http://www.madamedical.com)

MOR 6/05

## Specifications

---

<b>Supply Pressure:</b>	200-3,000 psi
<b>Outlet Pressure:</b>	High flow outlets: 45-55 psi at 200-3,000 psi inlet pressure
<b>Flow Capacity:</b>	High flow outlets: 100 lpm minimum single or combined flow with 200 psi inlet pressure Constant flow outlet selections: 0-25 lpm: 0, 1/4, 1/2, 1, 2, 3, 4, 6, 8, 10, 15, 25
<b>Filter:</b>	25 micron sintered bronze
<b>Relief Valve:</b>	Yes
<b>Operating Temperature:</b>	0° F to 125° F
<b>Storage Temperature:</b>	-30° F to 150° F
<b>Materials:</b>	Body: Nickel-plated brass Outlets: Nickel-plated brass Internal parts: Brass, stainless steel, Viton <sup>®</sup> , silicon and Teflon <sup>®</sup>

## Ordering Information

---

### **Yoke Pin Indexed CGA-870 Model Regulators For Small Oxygen Cylinders.**

<b>R1835GBHF-25</b>	1/4 -25 lpm with gauge, 2 check valves, hose barb outlet
<b>R1835-25BR</b>	1/2 -25 lpm with gauge, hose barb outlet
<b>1635-1</b>	Pressure reducer with gauge, 1 DISS male outlet (Not Shown)

### **CGA-540 Model Regulator For Large Oxygen Cylinders.**

<b>R1833-50-2</b>	Pressure reducer with gauge, 2 DISS check valve outlets
-------------------	---

**Additional models are available. For information about our complete line of oxygen regulators, contact your Mada customer service representative.**



### **MADA, INC.**

DOMESTIC/INTERNATIONAL  
625 Washington Avenue, Carlstadt, NJ 07072  
Tel: 201-460-0454 • Fax: 201-460-3509  
To Order Call Toll Free: **800-526-6370**  
Web Page: [www.madamedical.com](http://www.madamedical.com)

Distributed By: