

Mada Wall-Mounted Flowmeters

- *Precise administration in fractional increments.*
- *Available for oxygen or medical air.*
- *Adaptable for all 1/8 NPT inlet adapters.*
- *Back pressure compensated.*
- *Flush Flow.*
- *Inlet Filter.*



No. 1701 SERIES

No. 1701 Series:

- *Five year warranty.*
- *Durable cylindrical construction.*
- *Available in: 1/2 - 8 LPM, 1 - 15 LPM, 1/10 - 1 LPM*



MADA, INC.

DOMESTIC/INTERNATIONAL

625 Washington Avenue, Carlstadt, NJ 07072

Tel: 201-460-0454 • Fax: 201-460-3509

To Order Call Toll Free: **800-526-6370**

Web Site: www.madamedical.com

Mada Inlet Adapters

Mada provides a complete line of inlet adapters for the most popular oxygen and medical air flowmeters; factory assembled to your flowmeter.




No. 302-4C
Chrome Nut/Nipple



No. 302-4P
Plastic Nut/Nipple



No. 1703
Oxygen Fitting
DISS Female
(Nut & Nipple)
1/8" NPT

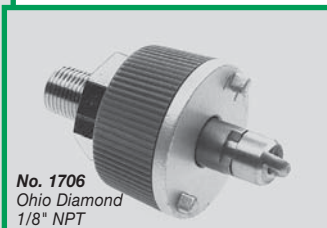


No. 1706
Ohio Diamond
1/8" NPT



No. 1708
Chemetron
(NCG) male 1/8" NPT



No. 1707
Puritan-Bennett
1/8" NPT

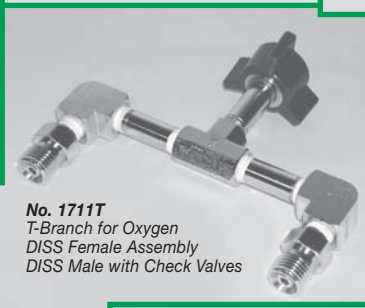


No. 1711T
T-Branch for Oxygen
DISS Female Assembly
DISS Male with Check Valves



No. 1714Y
Y-Branch for Oxygen
DISS Female Assembly
DISS Male with
Check Valves




No. 1720
Auxiliary Valve O₂ Power Take Off

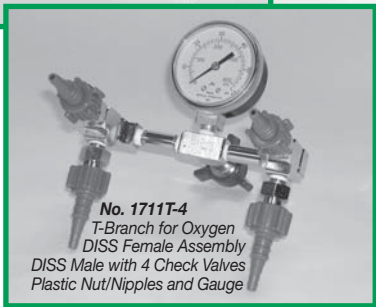

No. 1711T-4
T-Branch for Oxygen
DISS Female Assembly
DISS Male with 4 Check Valves
Plastic Nut/Nipples and Gauge

Adapters also available in medical air.

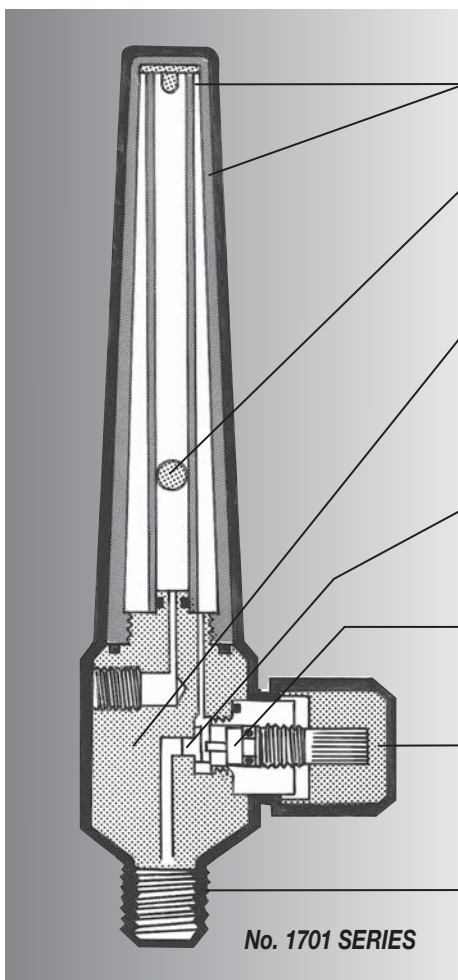

Dual scale calibrated flowtube and housing.

Polished, stainless steel float for excellent contrast when reading flow.

Chrome-plated, solid brass body for durability and beauty.

Metal seat for long life and stable flowrates.

Reliable, two-piece needle valve assembly.

Large ribbed knob for precision control.

DISS male outlet fitting.

No. 1701 SERIES

Specifications

MODEL	GAS SCALE	ACCURACY
1701	Oxygen 1-15 LPM	± 10%
1701-1	Oxygen 1/10-1 LPM	± 15%
1701-8	Oxygen 1/2-8 LPM	± 10%
1701-A	Air 1-15 LPM	± 10%

Distributed By:



MADA, INC.
DOMESTIC/INTERNATIONAL
625 Washington Ave., Carlstadt, NJ 07072
Tel: 201-460-0454 • Fax: 201-460-3509
Web Site: www.madamedical.com

To Order Call Toll Free: **800-526-6370**